

# Instructions for use

miniRITE  
miniRITE T

Oticon Jet



Made for  
iPhone | iPad | iPod

**oticon**  
life-changing technology



## Model overview

This booklet is valid for the Oticon Jet family in the following hearing aid models:

- miniRITE
- miniRITE T

### FW 1.0

- |                                       |                           |
|---------------------------------------|---------------------------|
| <input type="checkbox"/> Oticon Jet 1 | GTIN: (01) 05714464060962 |
| <input type="checkbox"/> Oticon Jet 2 | GTIN: (01) 05714464060979 |

The following speakers are available for the above models:

- Speaker 60
- Speaker 85
- Speaker 100 (Power Instrument)
  
- Power flex mould speaker 100 (Power Instrument)
- Power flex mould speaker 105 (Power Instrument)

## Introduction to this booklet

This booklet guides you on how to use and maintain your new hearing aid. Ensure you read this booklet carefully, including the **Warnings** section. This will help you get the most benefit from your hearing aid.

Your hearing care professional has adjusted the hearing aids to meet your needs. If you have additional questions, contact your hearing care professional.

A hearing care professional\* (hearing aid professional, audiologist, ENT (ear, nose and throat) doctor, and hearing aid dispenser) is a person who is appropriately educated and has proven competency in professionally assessing hearing, selecting, fitting, and delivering hearing instruments and rehabilitation care to persons with hearing loss.

\*The job title may vary from country to country.

The education of the hearing care professional is in accordance with national or regional regulations.

| **About** | [Getting started](#) | [Daily use](#) | [Options](#) | [Tinnitus](#) | [Warnings](#) | [More info](#) |

For your convenience, this booklet contains a navigation bar to help you easily navigate through the different sections.

## Intended use

<b>Intended use</b>	The hearing aid is intended to amplify and transmit sound to the ear.
<b>Indications for use</b>	Bilateral or unilateral impaired hearing of sensorineural, conductive, or mixed type ranging from a slight (16dB HL*) to profound (95dB HL*) degree of hearing loss, with an individual frequency configuration.
<b>Intended user</b>	Person with hearing loss using a hearing aid and their caregivers. Hearing care professional responsible for adjusting the hearing aid.
<b>Intended user group</b>	Adults and children older than 36 months.
<b>User environment</b>	Indoor and outdoor.
<b>Contraindications</b>	Not suitable for infants below 36 months. Users of active implants must pay special attention when using the hearing aid. For more information read the Warnings section.
<b>Clinical benefits</b>	The hearing aid is designed to provide better speech understanding to help ease communication with the aim of improving quality of life.

*\*As specified by the American Speech-Language-Hearing Association, [asha.org](http://asha.org), using pure-tone average of 0.5, 1 and 2 kHz.*

**IMPORTANT NOTICE**

The hearing aid amplification is uniquely adjusted and optimised to your personal hearing capabilities during the hearing aid fitting performed by your hearing care professional.

# Table of contents

## About

---

Your hearing aid, speaker, and earpiece	12
miniRITE	14
miniRITE T	15
Identify left and right hearing aid	16
MultiTool for handling batteries and cleaning	17

---

## Getting started

---

Turn hearing aid ON and OFF	18
When to replace the battery	19
How to replace the battery (size 312)	20
Put on the hearing aid	22

---

## Daily use

---

Caring for the hearing aid	23
Replace standard earpieces	24
Replace ProWax miniFit filter	27
Clean customised earpieces	28
Replace ProWax filter	29
Flight mode	30

---



## **Options**

Optional features and accessories	31
Change volume	32
Change program	34
Mute your hearing aids	36
Use hearing aid with iPhone and iPad	37
Wireless accessories	38
Other options	40

## **Tinnitus**

Tinnitus SoundSupport™ (optional)	42
Guidelines for Tinnitus SoundSupport users	43
Sound options and volume adjustment	44
Limitation on use time	48
Tinnitus SoundSupport warnings	51

## **Warnings**

General Warnings	52
------------------	----

Continues on next page


# Table of contents

## More info


Troubleshooting	58
Water & dust resistant (IP68)	60
Conditions of use	61
Technical information	62
Warranty	68
Your individual hearing aid settings	69
Sound indicators	71



# Your hearing aid, speaker, and earpiece

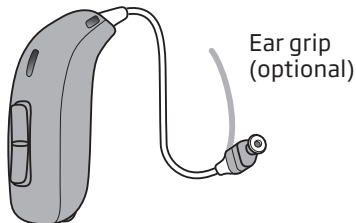

miniRITE T

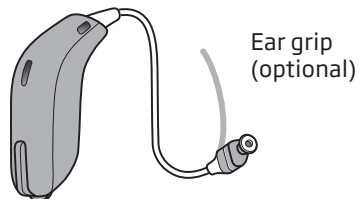


miniRITE


The hearing aid uses one of the following speakers:

**Standard speakers**


60

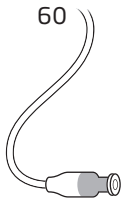

85

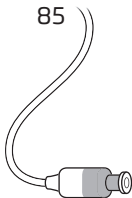


100

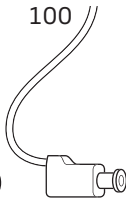

**Power flex mould speakers**

100

105

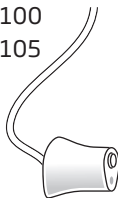

The speaker uses one of the following earpieces:

### Standard earpieces



Open dome



Bass dome, single vent



Bass dome, double vent



Power dome



Grip Tip

Available in small and large, left and right, with or without vent.

### Customised earpieces



Micro mould



LiteTip



VarioTherm® Micro mould



VarioTherm® LiteTip

® VarioTherm is a registered trademark of Dreve

### Dome sizes



5 mm\*

6 mm

8 mm

10 mm

12 mm\*\*


\* only as open dome for speaker 60

\*\* not for open dome

Please see details for replacing the dome in chapter "Replace standard earpieces".


# miniRITE

## What it is and does

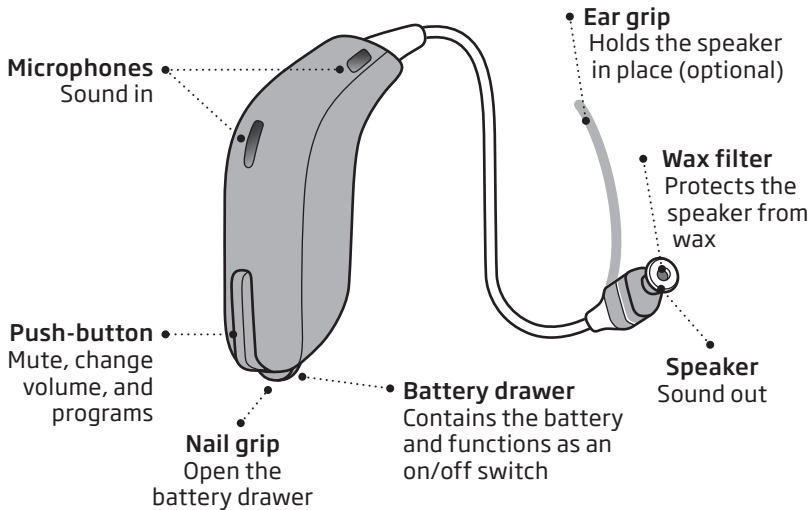

# miniRITE T

## What it is and does

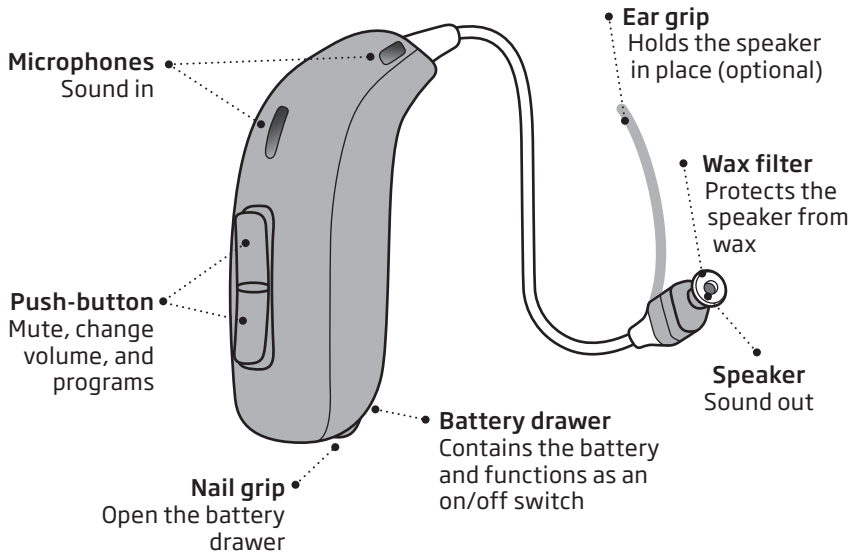


## Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid, as they may be programmed differently.

You can find left/right colour indicators in the battery drawer and on 60 and 85 speakers. Indicators (either L or R) can also be found on 100 speakers and some earpieces.


For 105 speakers, the indicator is found on the earpiece.

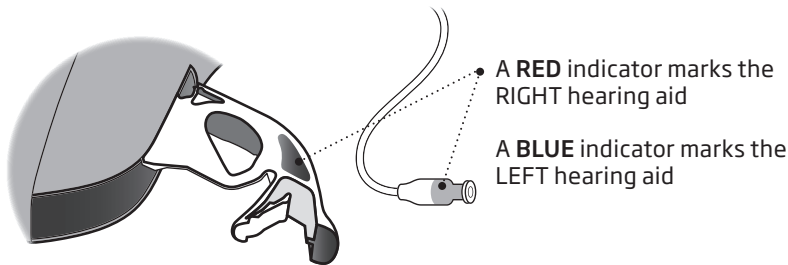



## MultiTool for handling batteries and cleaning

The MultiTool contains a magnet that makes it easier to replace the battery in the hearing aid. It also contains a brush and wire loop for cleaning and removing earwax. If you need a new MultiTool, please contact your hearing care professional.

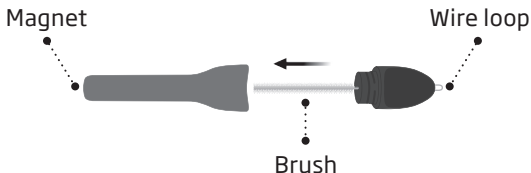

### IMPORTANT NOTICE


The MultiTool has a built-in magnet. Keep the MultiTool at least 30 centimetres (1 foot) away from credit cards and other magnetically sensitive devices.

## Turn hearing aid ON and OFF

The battery drawer is used to switch the hearing aid ON and OFF. To save battery life, make sure your hearing aid is switched OFF when you are not wearing it. If you wish to return to the standard settings of the hearing aid, simply open and close the battery drawer (quick reset).


### Turn ON

Close the battery drawer with the battery in place.

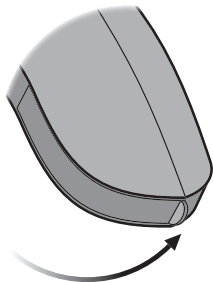


### Turn OFF


Open the battery drawer.

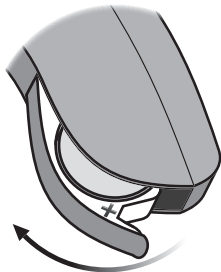

## When to replace the battery

When it is time to replace the battery, you will hear three tones repeated in moderate intervals until the battery runs out.

 **Three alternate tones\***  
= The battery is running low

 **Four descending tones**  
= The battery has run out

### Battery tip


To make sure the hearing aid is always working, bring spare batteries with you, or replace the battery before you leave home.

**Note:** Batteries need to be replaced more often if you are streaming to your hearing aid.

\* Bluetooth® will be turned off and it will not be possible to use wireless accessories.


# How to replace the battery (size 312)

## 1. Remove

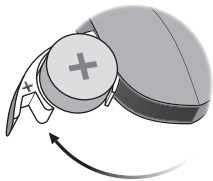

Fully open the battery drawer. Remove the battery.


## 2. Uncover

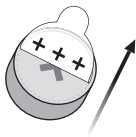

Remove the sticky label from the + side of the new battery.


Tip:  
Wait 2 minutes so that the battery can draw air, to ensure optimal functioning.

## 3. Insert

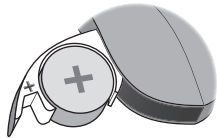

Insert the new battery into the battery drawer with the + side facing upwards.


#### 4. Close

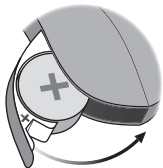

Close the battery drawer. The hearing aid will play a jingle through the earpiece.

Hold the earpiece close to your ear to hear the jingle.

#### Tip

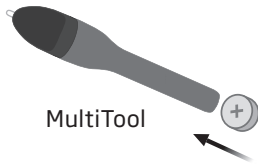

The MultiTool can be used to change the battery. Use the magnetic end to remove and insert batteries.

The MultiTool is provided by your hearing care professional.

# Put on the hearing aid


## Step 1



Place the hearing aid behind your ear.


The speaker should always be used with an earpiece attached. Use only parts designed for your hearing aid.

## Step 2

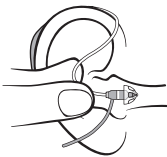

Hold the bend of the speaker wire between your thumb and index finger. The earpiece should point towards opening of the ear canal.

## Step 3

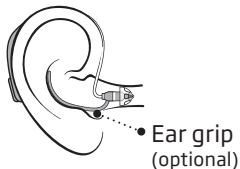

Gently push the earpiece into your ear canal until the speaker wire sits close to your head.

If the speaker has an ear grip, place it in the ear so it follows the contour of the ear.

## Caring for the hearing aid

When handling the hearing aid, hold it over a soft surface to avoid damage if you drop it.

### Clean the microphone openings

Use the brush of the MultiTool to carefully brush debris away from the openings. Carefully brush the surface around the opening. Make sure that no parts of the MultiTool are squeezed into the microphone openings by force. This may damage the hearing aid.




### IMPORTANT NOTICE

To clean the hearing aids, use a soft, dry cloth. The hearing aids must never be washed or immersed in water or other liquids.

## Replace standard earpieces


The standard earpiece (dome or Grip Tip) should not be cleaned. If the earpiece is filled with wax, replace it with a new one. Grip Tip should be replaced at least once a month.

### Step 1

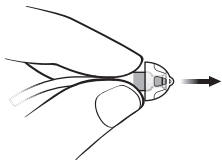


Hold on to the speaker and then pull off the earpiece.

### Step 2

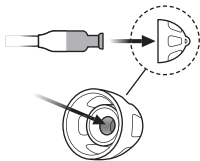

Insert the speaker exactly into the middle of the earpiece to obtain a secure attachment.

### Step 3

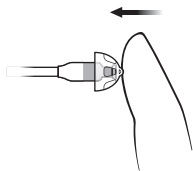

Push firmly to ensure that the earpiece is fastened securely.




**IMPORTANT NOTICE**

If the earpiece is not on the speaker when removed from the ear, the earpiece may still be in the ear canal. For further instructions, consult your hearing care professional.

## ProWax miniFit filter

The speaker has a white wax filter attached to the end where the earpiece is attached. The wax filter keeps wax and debris from damaging the speaker. Replace the filter when clogged, or if the hearing aid does not sound normal or contact your hearing care professional.

Remove the earpiece from the speaker before replacing the wax filter.

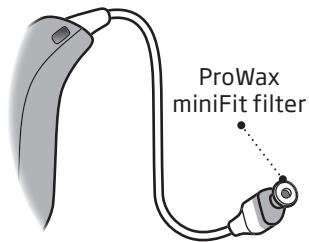


### IMPORTANT NOTICE

Ensure you always use the same type of wax filter as was originally supplied with the hearing aids. If you are in doubt about the use or replacement of wax filters, contact your hearing care professional.


# Replace ProWax miniFit filter

## 1. Tool

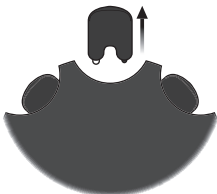


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

## 2. Remove

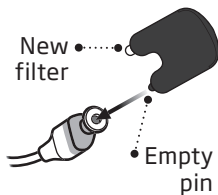

Push the empty pin into the wax filter in the speaker and pull it out.

## 3. Insert

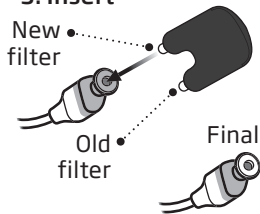

Insert the new wax filter using the other pin, remove the tool, and throw it out.

**Note:** If you use a mould or LiteTip, your hearing care professional must replace the wax filter in the speaker.


## Clean customised earpieces

The earpiece should be cleaned regularly.

The earpiece has a white wax filter.\* The filter keeps wax and debris from damaging the speaker.

Replace the filter when clogged, or if the hearing aid does not sound normal. Alternatively, contact your hearing care professional.


- Clean the vent by inserting the brush through the hole, twisting it slightly.

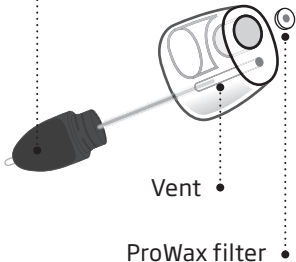

\* VarioTherm Micro mould and LiteTip do not have a wax filter


# Replace ProWax filter

## 1. Tool

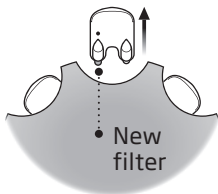


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

## 2. Remove

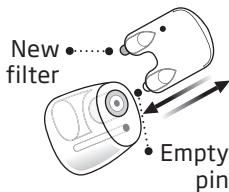

Push the empty pin into the wax filter in the earpiece and pull it out.

## 3. Insert

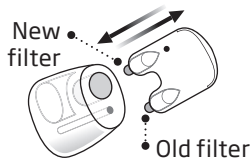


Insert the new wax filter using the other pin, remove the tool, and throw it out.


## Flight mode

When flight mode is activated, Bluetooth® is turned off. The hearing aid will still be working. Pressing the push-button on one hearing aid will activate flight mode on both hearing aids.

miniRITE

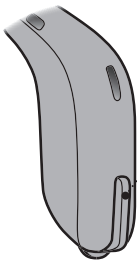

miniRITE T

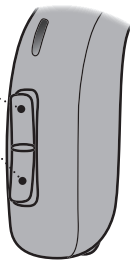

- **To activate and deactivate**  
Press the push-button for at least 7 seconds. A jingle confirms your action.

Opening and closing the battery drawer will also deactivate flight mode.

Press either end of the button

## Optional features and accessories

The features and accessories described in the following pages are optional. Please contact your hearing care professional for more information about features and accessories.

If you experience difficult listening situations, a special program may be helpful. These are programmed by your hearing care professional.

*Write down hearing situations in which you may need help.*

---


---

---

## Change volume

The push-button allows you to adjust the volume.  
You will hear a click when you increase or decrease the volume.

### miniRITE

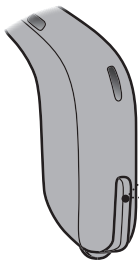


- A short press on the RIGHT hearing aid increases the volume
- A short press on the LEFT hearing aid decreases the volume





## miniRITE T

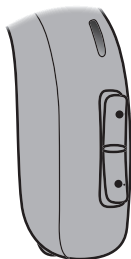

- Press button to increase volume

- Press button to decrease volume



## Change program


Your hearing aid can have up to four different programs. These are programmed by your hearing care professional. You will hear one to four tones when you change program depending on the program.

See the Sound indicator section.

## **miniRITE**

### **Press the button to change program**


If you have two hearing aids, the RIGHT hearing aid switches forwards, e.g. from program 1 to 2, and the LEFT hearing aid switches backwards, e.g. from program 4 to 3.

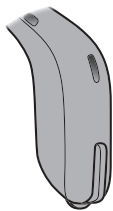

## **miniRITE T**

### **Press up or down between programs**


You can change continuously between programs - both up and down in the program order. For example, if you want to go from program 1 to 4, you can press the down button once instead of pressing the up button 3 times.

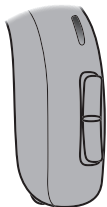

## Mute your hearing aids


Use the mute function if you need to silence the hearing aid. The mute function only mutes the microphone(s) on the hearing aid.

miniRITE

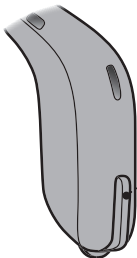

- Apply a long press (4 seconds) to the button to mute the hearing aid. To reactivate, push the button briefly.

miniRITE T

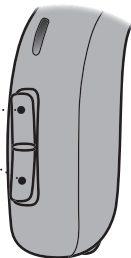

Press either end of the button

### IMPORTANT NOTICE

Do not use the mute function as an off switch, as the hearing aids are still using battery power in this mode.

## **Use hearing aid with iPhone and iPad**

Your hearing aids are Made for iPhone® and allow for direct communication and control with iPhone, iPad® or iPod touch®.

For assistance in using these products with your hearing aids, contact your hearing care professional or visit: [www.oticon.global/connectivity](http://www.oticon.global/connectivity)

To see how to pair your hearing aid with your iPhone, visit: [www.oticon.global/pairing](http://www.oticon.global/pairing)

For information on compatibility, visit: [www.oticon.global/compatibility](http://www.oticon.global/compatibility)

Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple products identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that use of this accessory with iPod, iPhone, or iPad may affect wireless performance.

## Wireless accessories

As an enhancement to your wireless hearing aid, a range of accessories is available. These can enable you to hear and communicate better in many everyday situations.

### ConnectClip

A device that can be used as remote microphone and handsfree headset when paired to your mobile phone.

### Phone Adapter 2.0

A device that when used together with hearing aids and ConnectClip, lets you communicate hands-free via a landline phone.

### TV Adapter 3.0

A device that streams sound from a TV or electronic audio device, to your hearing aids.

### Remote Control 3.0

A device that lets you change program, adjust volume, or mute your hearing aids.

## **Oticon ON**

An application that lets you control your hearing aid from your mobile phone or tablet. For iPhone, iPad, iPod touch, and Android devices. Ensure that you only download and install Oticon ON app from the official app stores.

For more information visit [www.oticon.global](http://www.oticon.global) or contact your hearing care professional.

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play, and the Google Play logo are trademarks of Google LLC.

## Other options



### **Telecoil - optional for miniRITE T**

Helps you hear better when using a telephone with a built-in loop, or when you are in buildings with teleloop systems such as theatres, churches, or lecture rooms. This symbol or a similar sign is shown wherever a teleloop has been installed.



### **Autophone - optional for Oticon Jet**

Activates a phone program in the hearing aid, if you place a dedicated magnet on your phone.

For more information, please contact your hearing care professional.





## **Tinnitus SoundSupport™ (optional)**

### **Intended use of Tinnitus SoundSupport**

Tinnitus SoundSupport™ is a tool intended to generate sounds to provide temporary relief for patients suffering from tinnitus as part of a tinnitus management program.

The target population is the adult population (over 18 years old).

Tinnitus SoundSupport is targeted to licensed hearing care professionals (audiologists, hearing aid specialists, or otolaryngologists) who are familiar with the evaluation and treatment of tinnitus and hearing loss. Fitting of Tinnitus SoundSupport must be done by a hearing care professional participating in a tinnitus management program.

## **Guidelines for Tinnitus SoundSupport users**

These instructions contain information about Tinnitus SoundSupport, which may have been enabled in your hearing aids by your hearing care professional.

Tinnitus SoundSupport is a tinnitus management device intended to generate sound of sufficient intensity and bandwidth to help manage tinnitus.

Your hearing care professional will also be able to offer the appropriate follow-up care. It is important to follow his/her advice and directions regarding such care.

### **Prescription use only**

Good health practice requires that a person reporting tinnitus has a medical evaluation by a licensed ear physician before using a sound generator. The purpose of such an evaluation is to ensure that any medically treatable condition that may cause tinnitus is identified and treated prior to using a sound generator.

## Sound options and volume adjustment

Tinnitus SoundSupport is programmed by your hearing care professional to match your hearing loss and preferences for tinnitus relief. It offers a number of different sound options. Together with your hearing care professional, you can select your preferred sound(s).

### Tinnitus SoundSupport programs

Together with your hearing care professional you decide for which programs you may want to have Tinnitus SoundSupport activated. The sound generator can be activated in up to four different programs.

### Mute

If you are in a program for which Tinnitus SoundSupport is activated, the mute functionality will mute only the environmental sounds, and not the sound from Tinnitus SoundSupport.

See chapter: “Mute your hearing aids”.

## **Volume adjustments with Tinnitus SoundSupport**

When you select a hearing aid program for which Tinnitus SoundSupport is activated, your hearing care professional can only set the push-button on your hearing aid to work as a volume control for the tinnitus relief sound.

Your hearing care professional will set the volume control for the sound generator in one of two ways:

- A) Change volume in each ear separately, or
- B) Change volume in both ears simultaneously.

## miniRITE

### **A) How to change Tinnitus SoundSupport volume in each ear separately**

To **increase** volume (on one hearing aid only), use a short press on the push-button repeatedly until desired level is reached. The sound will always be louder with the first press(es) until two beeps are heard. Hereafter the volume will decrease.

To **decrease** volume (on only one hearing aid), continue to press the push-button repeatedly until desired level is reached.

### **B) How to change Tinnitus SoundSupport volume in both ears simultaneously**

You can use one hearing aid to increase the sound and the other hearing aid to decrease the sound:

To **increase** volume, use a short press on the push-button repeatedly on the RIGHT hearing aid.

To **decrease** volume, use a short press on the push-button repeatedly on the LEFT hearing aid.

## miniRITE T

### **A) How to change Tinnitus SoundSupport volume in each ear separately**

To **increase** volume (on one hearing aid only), use a short press on the upper part of the push-button repeatedly until desired level is reached.

To **decrease** volume (on one hearing aid only), use a short press on the lower part of the push-button repeatedly until desired level is reached.

### **B) How to change Tinnitus SoundSupport volume in both ears simultaneously**

You can use one hearing aid to increase/decrease the sound in both hearing aids. When changing the volume in one hearing aid, the volume on the other hearing aid will follow.

To **increase** volume, use a short press on the upper part of the push-button repeatedly.

To **decrease** volume, use a short press on the lower part of the push-button repeatedly.

## Limitation on use time

### Daily use

The volume levels of Tinnitus SoundSupport can be set to a level which could lead to permanent hearing damage when used for a prolonged period of time. Your hearing care professional will advise you of the maximum amount of time per day you should use Tinnitus SoundSupport. It should never be used at uncomfortable levels.

See table “Tinnitus SoundSupport: Limitation on use” in section “Your individual hearing aid settings” at the end of this booklet to learn how many hours per day you can safely use the relief sound in your hearing aids.



# **Important information for hearing care professionals about Tinnitus SoundSupport**

## **Device description**

Tinnitus SoundSupport is a module function that can be enabled in the hearing aids by the hearing care professional.

## **Maximum wearing time**

The wearing time of Tinnitus SoundSupport will decrease as you increase the level above 80 dB(A) SPL. The fitting software will automatically display a warning when the hearing aid exceeds 80 dB(A) SPL. See “Max wearing time indicator” next to the tinnitus fitting graph in the fitting software.

## **The volume control is deactivated**

By default the volume control for the sound generator is deactivated in the hearing aid. Risk of noise exposure increases when the volume control is activated.

## **If the volume control is activated**

A warning may be displayed if you activate the tinnitus volume control in the “Buttons & Indicators” screen. This occurs if the relief sound can be listened to at levels that may cause hearing damage. The “Max wearing time” table in the fitting software displays the number of hours the patient can safely use Tinnitus SoundSupport.

- Note the max wearing time for each program for which Tinnitus SoundSupport is activated.
- Write those values in the table: “Tinnitus SoundSupport: Limitation on use”, in the back of this booklet.
- Instruct your patient accordingly.

## **Tinnitus SoundSupport warnings**

If your hearing care professional has activated the sound generator Tinnitus SoundSupport, please pay attention to the following warnings.

There are some potential concerns associated with the use of any sound generated by a tinnitus management device. Among them are the potential worsening of tinnitus, and/or a possible change in hearing thresholds.

Should you experience or notice a change in hearing or tinnitus, or any dizziness, nausea, headaches, heart palpitations, or possible skin irritation at the point of contact with the device, you should immediately discontinue use of the device and consult a medical, audiology, or other hearing care professional.

As with any device, misuse of the sound generator feature may cause potentially harmful effects. Care should be taken to prevent unauthorised use and to keep the device out of reach of children and pets.

### **Maximum wearing time**

Always follow the maximum wearing time per day of the Tinnitus SoundSupport advised by your hearing care professional. Prolonged use may lead to worsening of your tinnitus or of your hearing loss.

## **General warnings**

For your personal safety and to ensure correct usage, you should familiarise yourself fully with the following general warnings before using your hearing aids.

Consult your hearing care professional if you experience unexpected operations or serious incidents with your hearing aids during use or because of its use. Your hearing care professional will support you with issue handling and, if relevant, reporting to the manufacturer and/or the national authorities.

Note that hearing aids do not restore normal hearing and do not prevent or improve a hearing impairment resulting from organic conditions. Hearing aids are only a part of hearing habilitation and may need to be supplemented by auditory training and instruction in lipreading. Furthermore, note that in most cases, infrequent use of hearing aids does not permit a user to attain full benefit from it.

### **Usage of hearing aids**

Hearing aids should be used only as directed and adjusted by your hearing care professional. Misuse can result in sudden and permanent hearing loss.

Never allow others to wear your hearing aid as incorrect usage could cause permanent damage to their hearing.

### **Choking hazards & risk of swallowing batteries and other small parts**

Hearing aids, their parts and batteries should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves.

Batteries have occasionally been mistaken for pills. Therefore, check your medicine carefully before swallowing any pills.

*If a battery, hearing aid or small part is swallowed, see a doctor immediately.*

### **Battery use**

Always use batteries recommended by your hearing care professional. Batteries of low quality may leak and cause bodily harm.

Never attempt to recharge your batteries and never dispose of batteries by burning them. There is a risk that the batteries will explode.

### **Dysfunction**

Be aware of the possibility that your hearing aids may stop working without notice. Keep this in mind when you depend on warning sounds (e.g. when you are in traffic). The hearing aids may stop functioning for instance if the batteries have expired or if the tubing is blocked by moisture or earwax.

### **Active implants**

The hearing aid has been thoroughly tested and characterised for human health according to international standards for human exposure (Specific Absorption Ratio - SAR), induced electromagnetic power and voltages into the human body.

The exposure values are well below international accepted safety limits for SAR, induced electromagnetic power and voltages into the human body defined in the standards for human health and coexistence with active medical implants, such as pacemakers and heart defibrillators.

If you have an active brain implant, contact the manufacturer of your implantable device for information about the risk of disturbance.

## **General warnings**

The Autophone magnet or MultiTool (which has a built-in magnet) should be kept more than 30 centimetres (1 foot) away from the implant, e.g. do not carry it in your breast pocket. Follow the guidelines recommended by the manufacturers of implantable defibrillators and pacemakers regarding their use with magnets.

### **Cochlear implants**

If you are using a cochlear implant (CI) on one ear and a hearing aid on the other ear, make sure to always keep your CI at least at a 1 centimetre (0.4 inches) distance from your hearing aid. The magnetic field from the CI sound processors, coils, and magnets may permanently damage the speaker unit in your hearing aid. Never place the devices close together on a table e.g. when cleaning or changing batteries. Do not carry the CI and the hearing aid together in the same box.

### **X-ray/CT/MR/PET scanning, electrotherapy and surgery**

Remove your hearing aid before X-ray examinations and CT/MR/PET scans, electrotherapy, surgery, etc. as your hearing aid may be damaged when exposed to strong electromagnetic fields.

### **Heat and chemicals**

The hearing aid must never be exposed to extreme heat, e.g. left inside a parked car in the sun.

The hearing aid must not be dried in microwave ovens or other ovens.

The chemicals in cosmetics, hairspray, perfume, aftershave lotion, suntan lotion, and insect repellent can damage the hearing aid. Always remove your hearing aid before applying such products and allow time to dry before use.

**Detached earpiece in ear canal**

If the earpiece is not on the speaker when removed from the ear, the earpiece may still be in the ear canal. For further instructions, consult your hearing care professional.

**Power hearing aid**

Special care should be exercised in selecting, fitting, and using hearing aids where the maximum sound pressure capability exceeds 132 dB SPL (IEC 60318-4 / IEC 711) as there may be a risk of impairing the remaining hearing of the hearing aid user.

*For information on whether your instrument is a power hearing aid, see the model overview section in this booklet.*

**Possible side effects**

Hearing aids, moulds or domes may cause an accelerated accumulation of earwax.

The non-allergenic materials used in hearing aids may, in rare cases cause a skin irritation or other side effects.

*If these conditions occur, seek consultation with a physician.*

**Use on aircraft**

Your hearing aids have Bluetooth wireless technology. On board an aircraft, the hearing aids must be put into Flight mode to deactivate Bluetooth, unless Bluetooth is permitted by the flight personnel.

## **General warnings**

### **Connection to external equipment**

The safety of the hearing aids, when connected to external equipment (with an auxiliary input cable and/or USB cable and/or directly), is determined by the external signal source. When the device is connected to external equipment which is plugged into a power socket, this equipment must comply with IEC 62368-1 or equivalent safety standards.

### **Use of third-party accessories**

Only use accessories, transducers or cables supplied by the manufacturer. Non-original accessories may result in reduced electromagnetic compatibility (EMC) of your device.

### **Modification of hearing aids is not allowed**

Changes or modifications not expressly approved by the manufacturer will void the warranty of the equipment.

### **Interference**

The hearing aids have been thoroughly tested for interference according to the most stringent international standards.

Electromagnetic interference may occur in the vicinity of equipment with the symbol to the left. Portable and mobile RF (radio frequency) communications equipment can affect the performance of your hearing aids. If your hearing aids are affected by electromagnetic interference, move away from the source to reduce the interference.





# Troubleshooting

Symptom	Possible causes
<b>No sound</b>	Dead battery
	Clogged earpieces (dome, Grip Tip, or mould)
	Hearing aid microphone muted
<b>Intermittent or reduced sound</b>	Clogged sound outlet
	Moisture
<b>Squealing noise</b>	Hearing aid earpiece inserted incorrectly
	Earwax accumulated in ear canal
<b>Beeping</b>	If your hearing aid plays 8 beeps, 4 times consecutively, your hearing aid needs a microphone service check
<b>Pairing issue with Apple device</b>	Bluetooth connection failed
	Only one hearing aid paired

## Solutions

Replace the battery

Clean mould

Replace wax filter, dome, or Grip Tip

Unmute the hearing aid microphone

Clean mould or replace wax filter, dome, or Grip Tip

Wipe battery with a dry cloth

Re-insert the earpiece

Have ear canal examined by your doctor

Contact your hearing care professional

- 1) Unpair your hearing aid.
- 2) Turn Bluetooth off and on again.
- 3) Open and close battery drawer on hearing aid.
- 4) Re-pair hearing aid (please visit [www.oticon.global/pairing](http://www.oticon.global/pairing)).

If none of the above solutions work, consult your hearing care professional.

## Water & dust resistant (IP68)

Your hearing aid is dust tight and protected against ingress of water, which means it is designed to be worn in all daily life situations. Therefore, you do not have to worry about sweat or getting wet in the rain. Should your hearing aid come into contact with water and stop working, please follow these guidelines:

1. Gently wipe off any water.
2. Open the battery drawer and remove the battery and gently wipe off any water in the battery drawer.
3. Let the hearing aid dry with the battery drawer left open for approximately 30 minutes.
4. Insert a new battery.

### IMPORTANT NOTICE

Do not wear your hearing aid while showering or participating in water activities. Do not immerse your hearing aid in water or other liquids.

## Conditions of use

<b>Operating conditions</b>	Temperature: +1°C to +40°C (34°F to 104°F) Humidity: 5% - 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
<b>Transportation and storage conditions</b>	Temperature and humidity shall not exceed the below limits for extended periods during transportation and storage:  <b>Transportation:</b> Temperature: -25°C to + 60°C (-13°F to 140°F) Humidity: 5% - 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa <b>Storage:</b> Temperature: -25°C to + 60°C (-13°F to 140°F) Humidity: 5% - 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa

## Technical information

The hearing aids contain the following two radio technologies:

The hearing aids contain a radio transceiver using short range magnetic induction technology operating at 3.84 MHz. The magnetic field strength of the transmitter is very weak and always below 15 nW (typically below  $-40 \text{ dB}\mu\text{A}/\text{m}$  at a distance of 10 meters ( $-12.20 \text{ dB}\mu\text{A}/\text{ft}$ ) at a distance of 33 feet).

The hearing aids also contain a radio transceiver using Bluetooth Low Energy technology and a proprietary, short-range radio technology both operating at ISM band 2.4 GHz.

The radio transmitter is weak and always below 3 mW equal to 4.8 dBm in total radiated power.

Only use your hearing aids in areas where wireless transmission is permitted.

The hearing aids comply with international standards concerning radio transmitters, electromagnetic compatibility, and human exposure.

Due to the limited space available on the hearing aid, relevant approval markings can be found in this booklet.

Additional information can be found in the "Technical Data Sheet" on [www.oticon.global](http://www.oticon.global)

The manufacturer declares that this hearing aid is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

This medical device complies with Medical Device Regulation (EU) 2017/745.

Declaration of Conformity is available at the headquarters.

Oticon A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark  
[www.oticon.global/doc](http://www.oticon.global/doc)

CE 0123



SBO Hearing A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark



Waste from electronic equipment must be handled according to local regulations.



 **Bluetooth®**

**IP68**

## Description of symbols and abbreviations used in this booklet



### Warnings

Text marked with a warning symbol must be read before using the device.



### Manufacturer

The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Regulation 2017/745.



### CE mark

The device complies with all required EU regulations and directives. The four digit number indicates the identification of the notified body.



### Electronic waste (WEEE)

Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).



### Regulatory Compliance Mark (RCM)

The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand markets.

IP68

### IP code

Class of protection against harmful ingress of water and particulate matter according to EN 60529. IP6X indicates total dust protection. IPX8 indicates the protection against the effects of continuous immersion in water.



### Bluetooth logo

Registered trademark of Bluetooth SIG, Inc. where any use of such requires a license.



Available on  
iPhone | iPad | iPod

### **Made for Apple badges**

The device is compatible with iPhone, iPad and iPod touch.



### **Hearing loop**

This logo incorporates the universal symbol for hearing assistance. The "T" signifies that a hearing loop is installed.



### **Radio Frequency (RF) transmitter**

Your hearing aid contains an RF transmitter.

**GTIN**

### **Global Trade Item Number**

A globally unique 14-digit number used to identify medical device products including medical device software. GTIN in this booklet is related to medical device firmware (FW). GTIN on regulatory packaging label is related to medical device hardware.

**FW**

### **Firmware**

Firmware version used in the device.

## Description of additional symbols used on labels



### **Keep dry**

Indicates a medical device that needs to be protected from moisture.



### **Catalogue number**

Indicates the manufacturer's catalogue number so that the medical device can be identified.



### **Serial number**

Indicates the manufacturer's serial number so that a specific medical device can be identified.



### **Medical Device**

The device is a medical device.



### **Unique device identifier**

Indicates a carrier that contains unique device identifier information.

## **International warranty**

Your device is covered by an international warranty issued by the manufacturer. This international warranty covers manufacturing and material defects in the device itself, but not in accessories such as batteries, tubing, speakers, earpieces and filters, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, physical changes in your ear, damage due to foreign objects entering the device, or incorrect adjustments are NOT covered by the international warranty and may void it. The above international warranty does not affect any legal rights that you might have under applicable national legislation governing the sale of consumer goods in the country where

you have bought your device. Your hearing care professional may also have issued a warranty that goes beyond the clauses of this international warranty. Please consult him/her for further information.

### **If you need service**

Take your device to your hearing care professional, who may be able to sort out minor problems and adjustments immediately. Your hearing care professional may charge a fee for their services.

# Warranty

## Certificate

Name of owner: \_\_\_\_\_

Hearing care professional: \_\_\_\_\_

Hearing care professional's address: \_\_\_\_\_

Hearing care professional's phone: \_\_\_\_\_

Purchase date: \_\_\_\_\_

Warranty period: \_\_\_\_\_ Month: \_\_\_\_\_

Model left: \_\_\_\_\_ Serial no.: \_\_\_\_\_

Model right: \_\_\_\_\_ Serial no.: \_\_\_\_\_

# Your individual hearing aid settings

*To be filled out by your hearing care professional.*

Tinnitus SoundSupport: Limitation on use			
<input type="checkbox"/>	No limitation on use		
	Program	Start-up volume (Tinnitus)	Max volume (Tinnitus)
<input type="checkbox"/>	1	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	2	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	3	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	4	Max _____ hours per day	Max _____ hours per day

Continues on next page

## Settings overview for your hearing aid

Left			Right	
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Change volume	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Change program	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Mute	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Tinnitus SoundSupport	<input type="checkbox"/> Yes	<input type="checkbox"/> No
Volume control indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at min/max volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps when changing volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at start-up volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
Battery indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Low battery warning	<input type="checkbox"/> On	<input type="checkbox"/> Off

## Sound indicators

Different sounds indicate the hearing aid status.  
The different indicators are listed on the following pages.

Your hearing care professional can set sound indicators to match your preferences.

Program	<input type="checkbox"/> Sound	When to use
1	1 tone	
2	2 tones	
3	3 tones	
4	4 tones	

Continues on next page

<b>ON</b>	<b>Sound</b>
On	<input type="checkbox"/> jingle
<b>Volume</b>	<b>Sound</b>
Start-up volume	<input type="checkbox"/> 2 beeps
Minimum/maximum volume	<input type="checkbox"/> 3 beeps
Volume up/down	<input type="checkbox"/> 1 beep
Mute	

<b>Accessories</b>	<input type="checkbox"/> <b>Sound</b>
TV Adapter	2 different tones
ConnectClip remote microphone	2 different tones
<b>Flight mode</b>	<b>Sound</b>
Flight mode active	Short jingle
Flight mode inactive	Short jingle

<b>Warnings</b>	<b>Sound</b>
Low battery	<input type="checkbox"/> 3 alternate tones
Battery shut down	4 descending tones
Microphone service check needed	8 beeps repeated 4 times







**oticon**  
life-changing **technology**